

GMFCS
Levels I-V


I am your faithful recumbent trike that grows with you.

Let's be friends for life! Friends always support each other and are there when it matters – and that's exactly what I'll do for you. With my growing system, I can accompany you from a young age and stay by your side as you grow. And because moving is so important – and feels so good, I can't wait for all the adventures we'll share. You'll be surprised how smoothly you can pedal with me. My special seat helps you bend your hips and knees comfortably thanks to a stable sitting position, so you can reach your goals with little effort. Hop on, and let's discover the world together!



Adjustable wheelbase


enables a long
Service life for users with
a height of 110–175 cm.

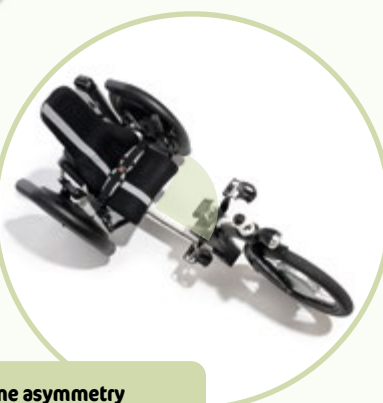

Frame asymmetry

provides free access to the
transfer vehicle and enables
easy and independent
boarding and alighting.



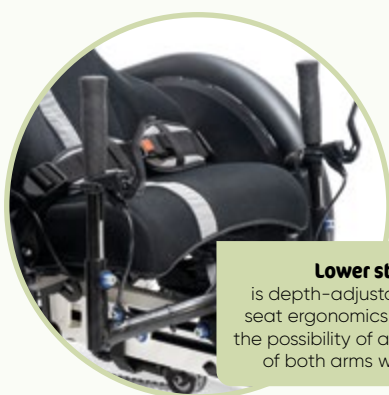

Lower steering
is depth-adjustable to suit the seat ergonomics and thus offers the possibility of a relaxed position of both arms when steering.



Rain cape
protects the seat from the weather.

Frame colours

-  Pure white
-  Turquoise-blue
-  Pure orange



Let's go on an adventure and discover the world together! I will accompany you as you grow with my ingenious and growing frame.



Automatic chain tensioner
effectively prevents the chain
from sagging.



Ergonomic seat
also offers a very low seating
position and, in combination
with the special frame geometry,
ensures optimum cornering.

data	mats.
Height	110–175 cm
Seat bone length	approx. 65–110 cm
Inside leg length	46–83 cm
Bottom bracket height	32.5 cm
Entry level	20 cm
Seat height	32–42 cm
Seat angle	0–20°
Seat depth	26 cm
Length of the backrest	50 cm
Wheelbase	115–142 cm
Track width	69 cm
Total length	165–192 cm
Total width	76 cm
Total height	65–88 cm
Max. seat width narrow seat	37 cm
Max. seat width wide seat	46 cm
Turning circle	243–395 cm
Tare weight	29 kg
Max. load	75 kg
HMV-No.	22.51.02.0063