

GMFCS
Level III-V

kg
max. 150



cm
26–46

cm
31–48

°
(-15° – 30°)



I am your comfortable shower and toilet chair.

Whether cleaning your teeth or showering, my seat tilt can be adapted to any situation. My anatomically shaped seat surface gives you the required support to stabilise your upper body and can be adjusted in height. You can sit comfortably on my soft PU pads which can be easily removed and cleaned. Because I have been produced from highly robust material, you and your parents will have a lot of fun with me for a long time to come. Did you know that I am available in five sizes and can carry up to 150 kg?



You can sit really comfortably on the **removable seat pad**.



Swing away arm supports also help to stabilise you.



Headrest
can be adjusted in
height and angle.

Tray
with padded armrest.

Insertion pad
prevents pressure marks on the body.

Seat pan
Anatomically moulded provides
stability in the pelvic area and
stabilises the upper body.

Robust sub-frame
made from stainless steel.

Footplate
can simply be folded up
for easier access.


Foot pedal
adjusts the seat pan continuously
from (-)5° to 30°.



Toilet equipment

You can use me as a toilet
chair independently of a WC.

Sizes 2-4
also feature a
larger hygiene
opening.

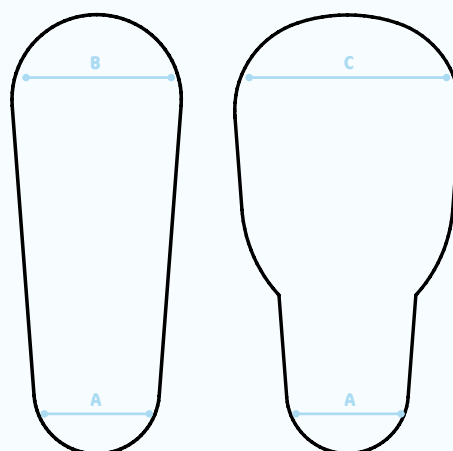



data		Size 0	Size 1	Size 2	Size 3	Size 4
Seat width (rear)		26 cm	31 cm	36 cm	41 cm	46 cm
Seat depth		31 cm	36 cm	41 cm	46 cm	48 cm
Seat shell width (front)		31 cm	36 cm	41 cm	46 cm	51 cm
Back height		49 cm	51 cm	56 cm	61 cm	65 cm
Seat angle		(-)5° - 30°				
Can be positioned over a toilet		47 - 57 cm	45 - 55 cm	42 - 52 cm		
Total width		56 cm	61 cm	63 cm	71,5cm	74 cm
Total length		65 cm			70 cm	
Max. lower leg length		38 cm		48 cm		
Front edge of hygiene opening to rear edge of base frame		47 cm	49 cm	51 cm	54 cm	57 cm
Wide opening at the front	A	10 cm				
Wide opening at the rear	B	12 cm		13,5cm		
Wide opening at the rear*	C	-	-	18,5cm		
Max. load		60 kg	80 kg	100 kg		150 kg
Weight		28 kg	30 kg	33 kg	35 kg	38 kg
HMV-No.		33.40.04.1014				

* Applies only to sizes 2b, 3b and 4b



You can find more information on accessories & harnesses in the price list.

